


NGS - quality control, alignment, visualisation


Sequencing technologies


Major applications

- Transcriptome characterization
 - e.g. RNA-seq
- Epigenome characterization:
 - e.g. ATAC-seq
- DNA-protein interactions:
 - e.g. ChIP-seq
- Whole genome (assembly)
- Variant detection
- Metagenome characterization
- Any others?

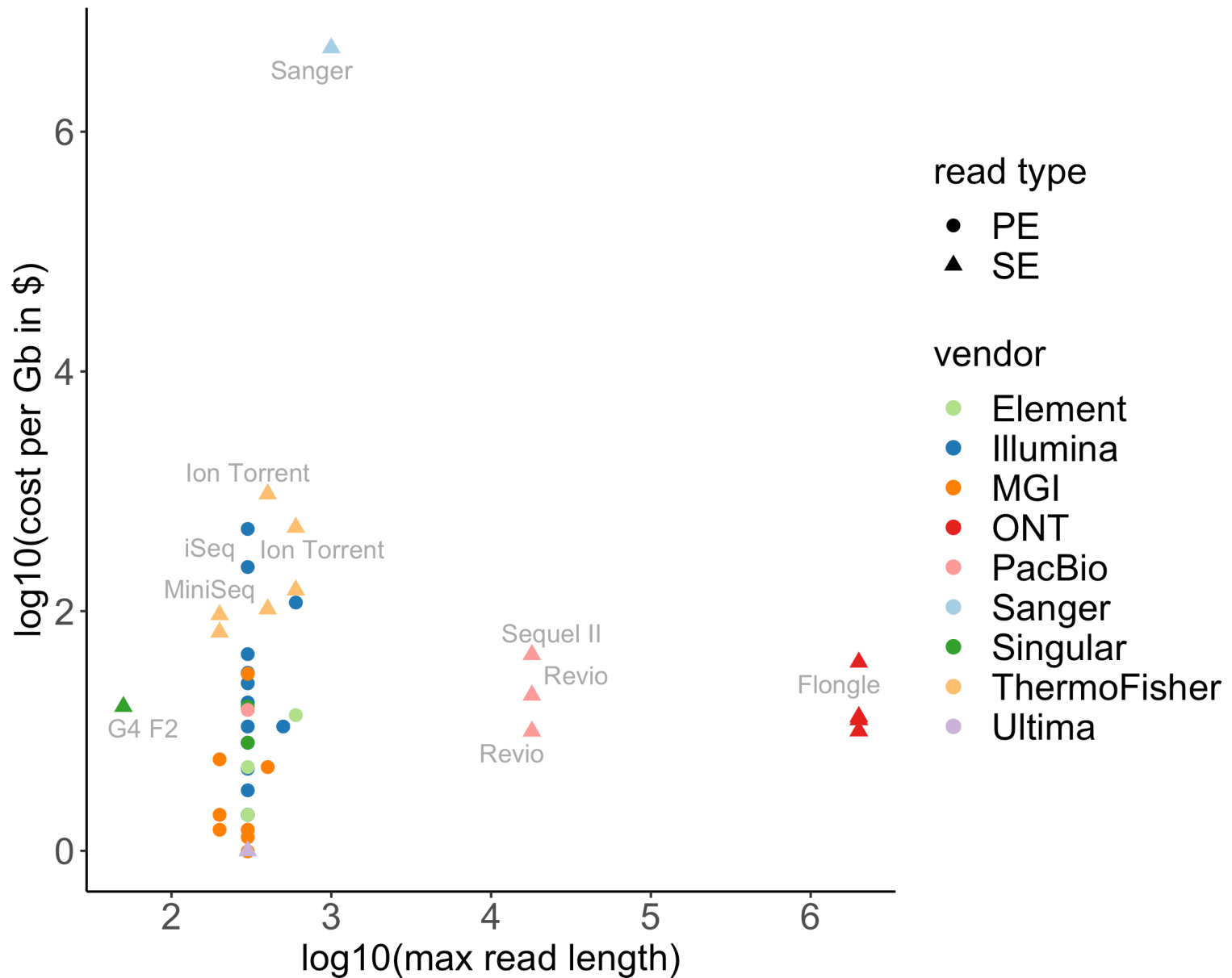

Quiz Question 4


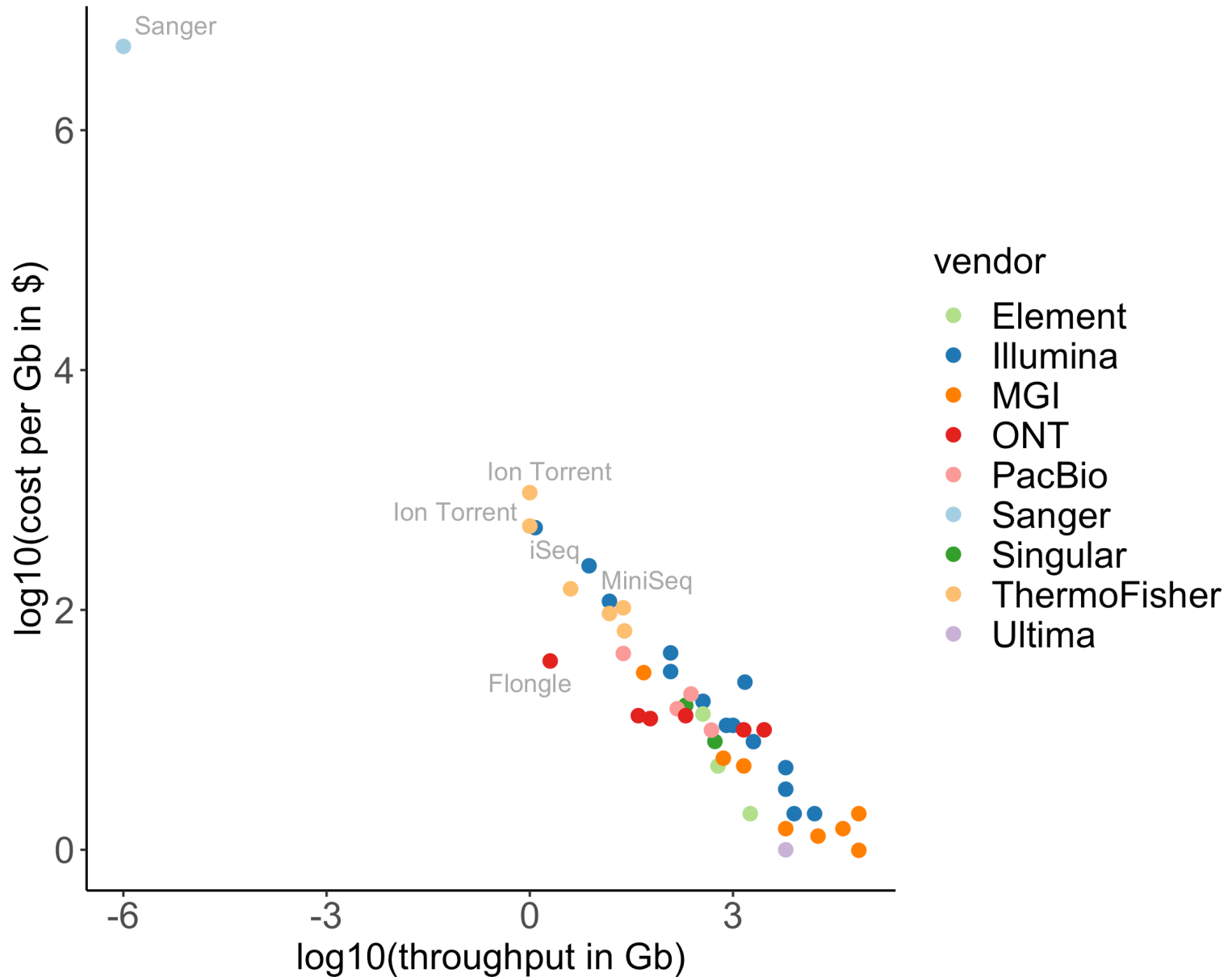


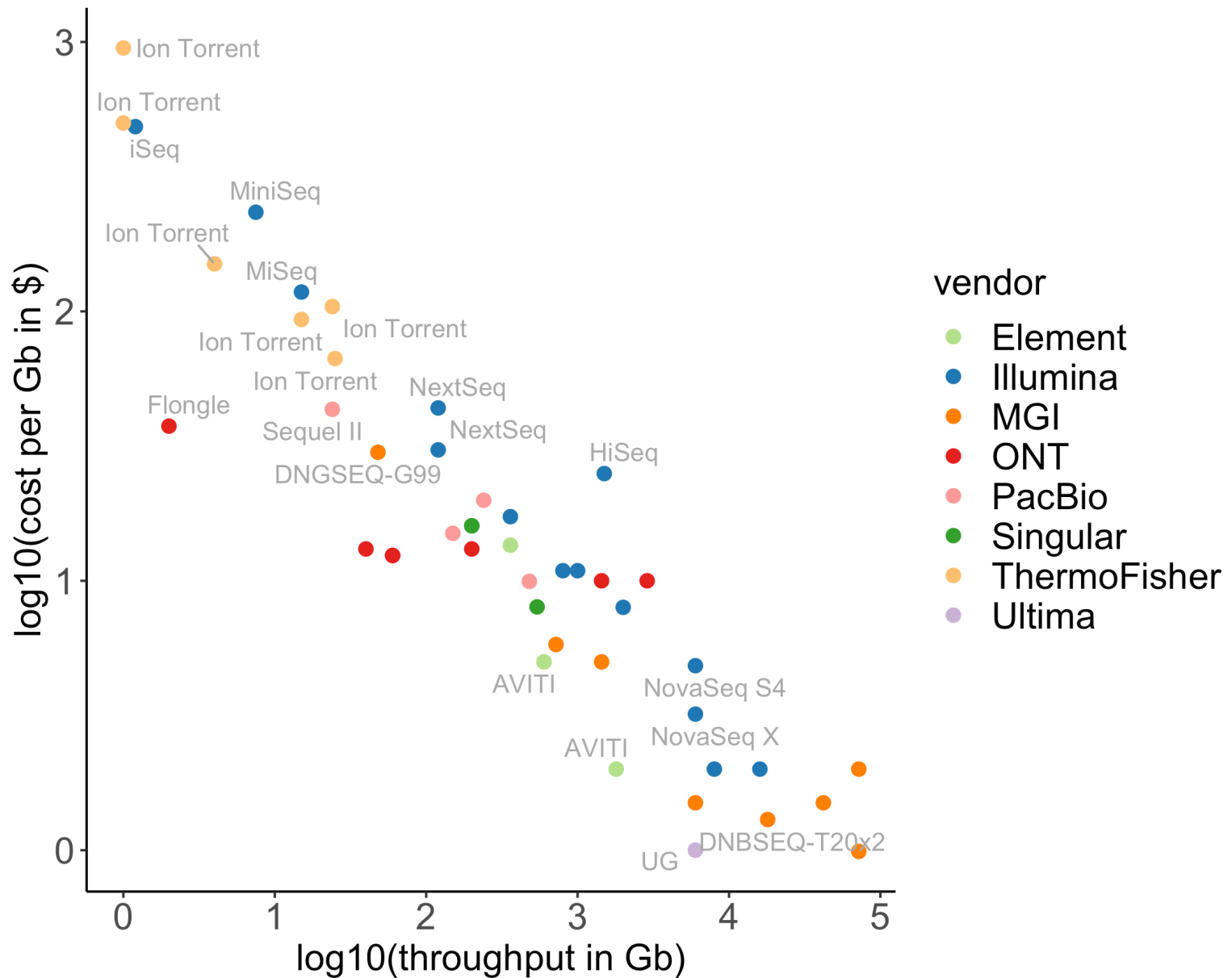

This course

- 2nd generation (sequencing by synthesis):
 - Ion Torrent
 - Illumina
- 3rd generation:
 - Pacific Biosciences
 - Oxford Nanopore Technology

Ion Torrent sequencing

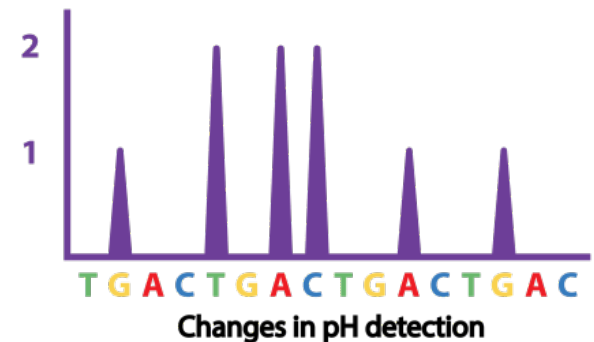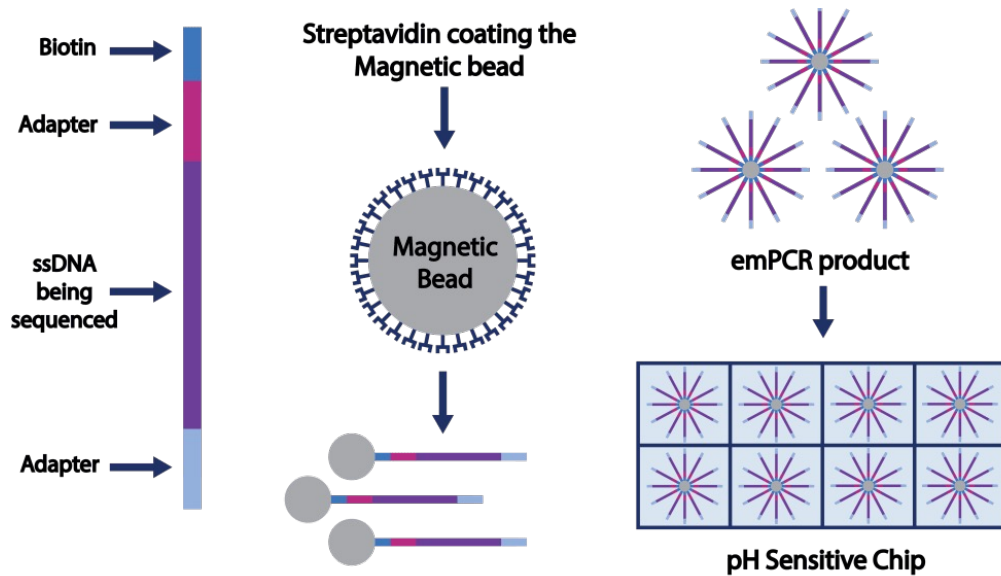

Ion Torrent sequencing

- Up to ± 400 bp read length
- Scalable (but Illumina has similar size systems nowadays)
- Homopolymers (e.g. TTTT) are a challenge (impossible) to sequence




Illumina sequencing

- Massive throughput: up to 16×10^{12} bases/run (NovaSeq X)
- Most used platform today

Illumina sequencing

- 50 – 300 bp
- Paired-end (or single-end)
- Multiplexing

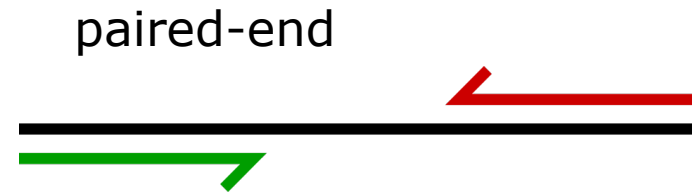

Illumina library prep

shear + size select DNA



Ligate adapters




Barcode + p5/p7 sites




PCR: 8-16 cycles




Sequencing

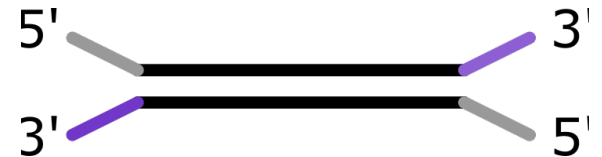

Sequencing by synthesis

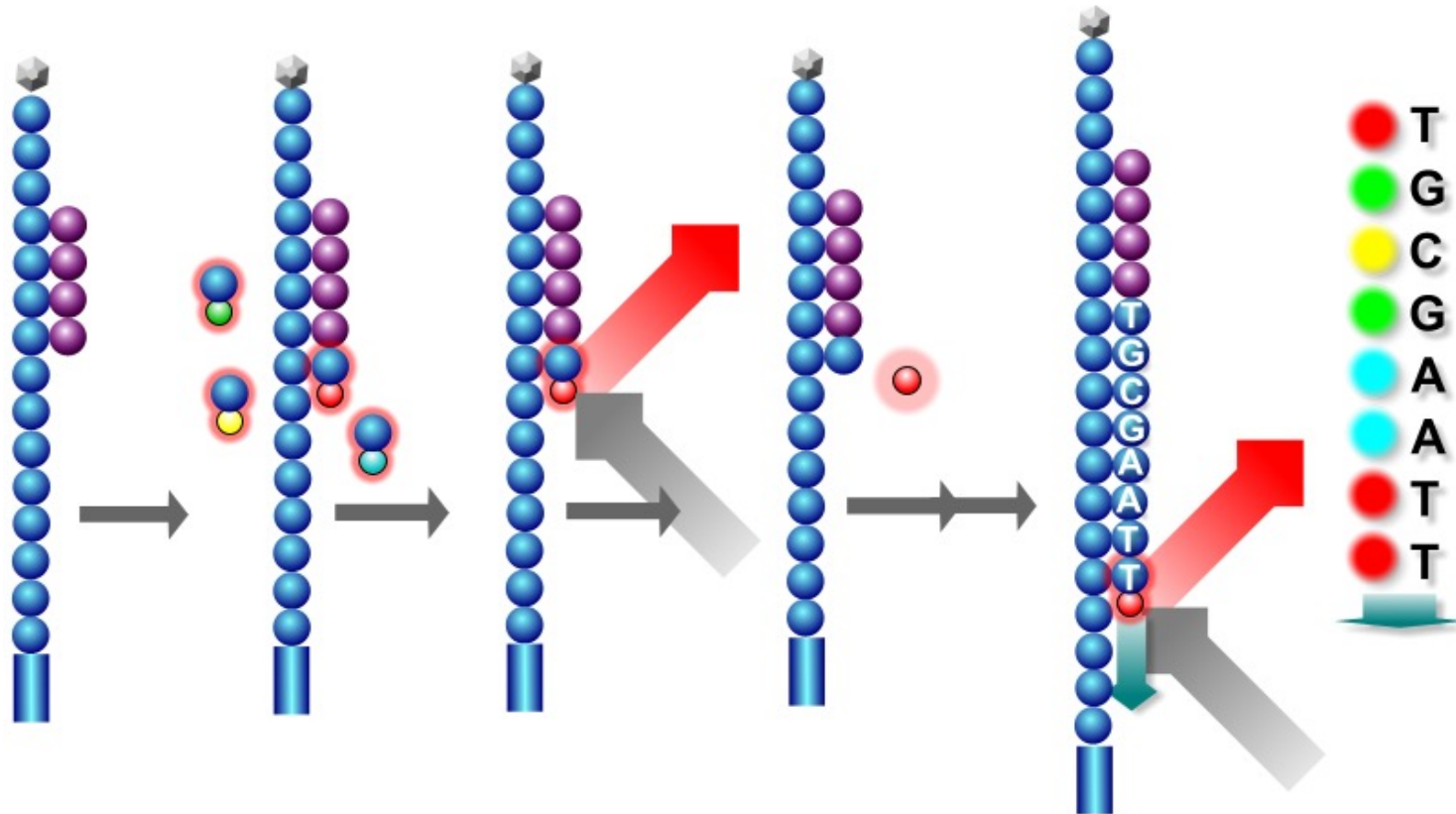

Some definitions

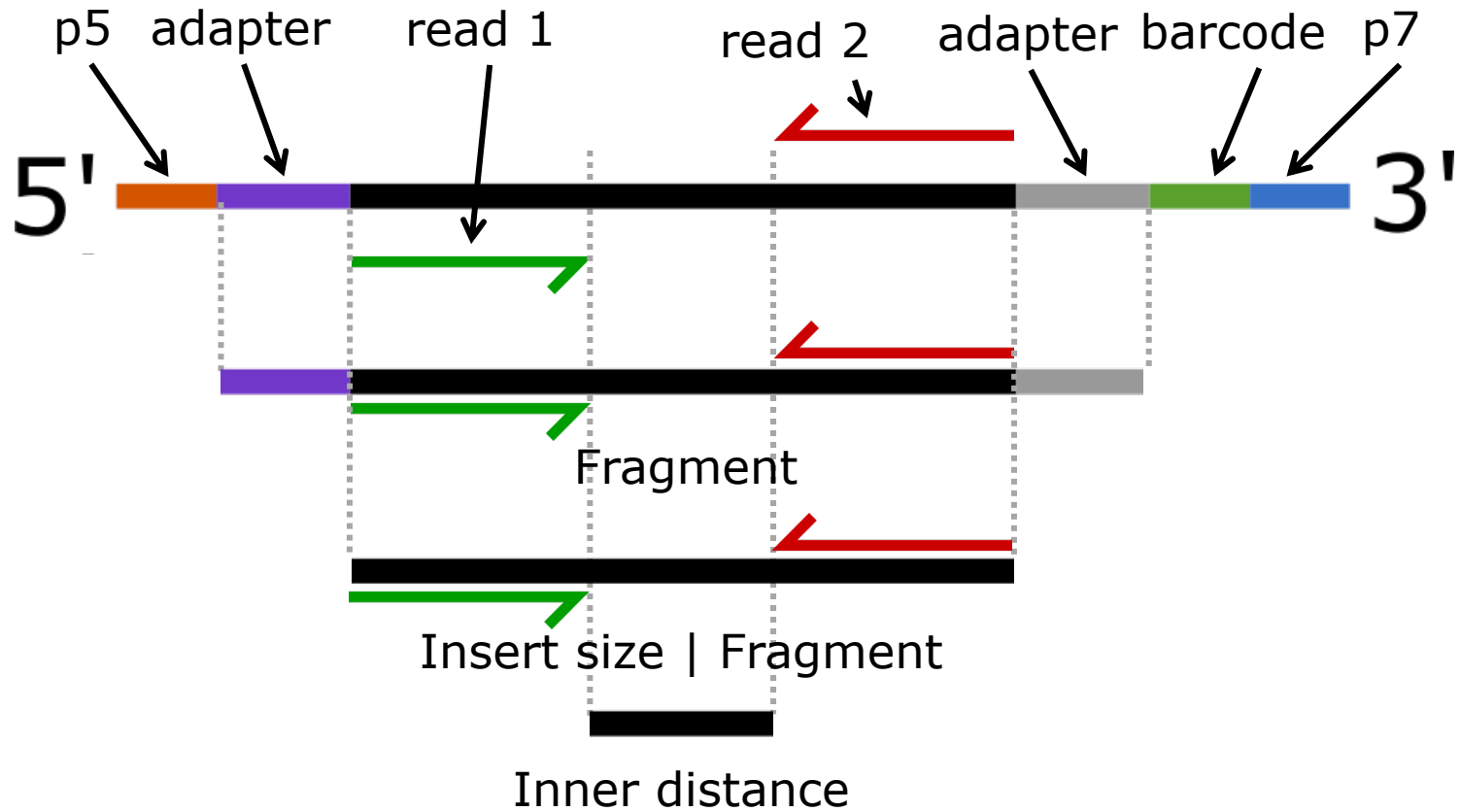

Some more definitions..

- **Library:** fragments from one (c)DNA sample that share a barcode
- **Sequencing run:** complete cycle of generating reads on a machine
- **Flow cell:** physical platform where sequencing reactions take place. Used once in a sequencing run.
- **Lane:** compartment within the flow cell. An Illumina flow cell often has multiple lanes (2 or 4)


Quiz Question 5



Illumina - limitations

- Maximum read length: 300 bp
- How to reconstruct:
 - Repeats?
 - Isoforms?
 - Structural variation?
 - Haplotypes?
 - Genomes?
- Why not longer read lengths?


Illumina - limitations

Sequence quality declines towards the end

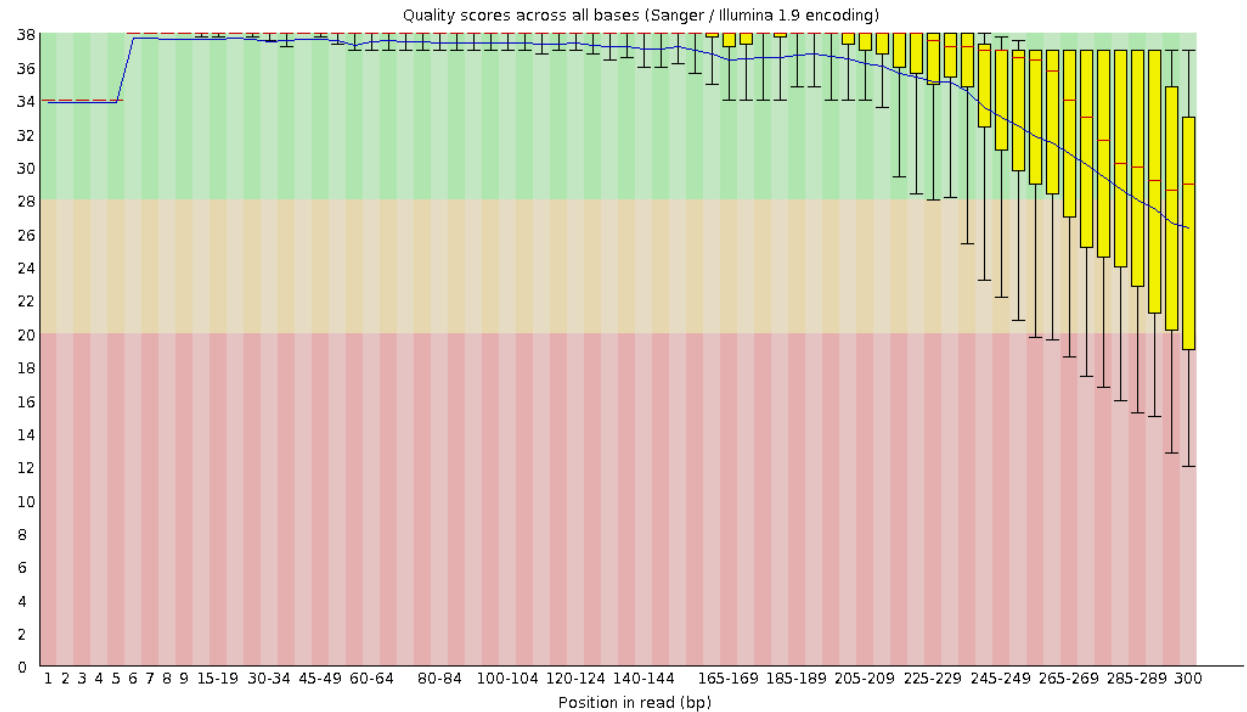

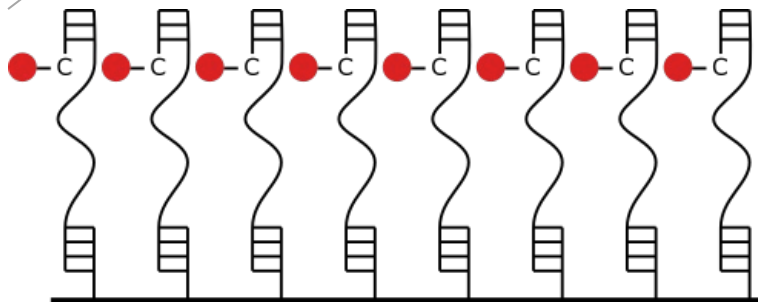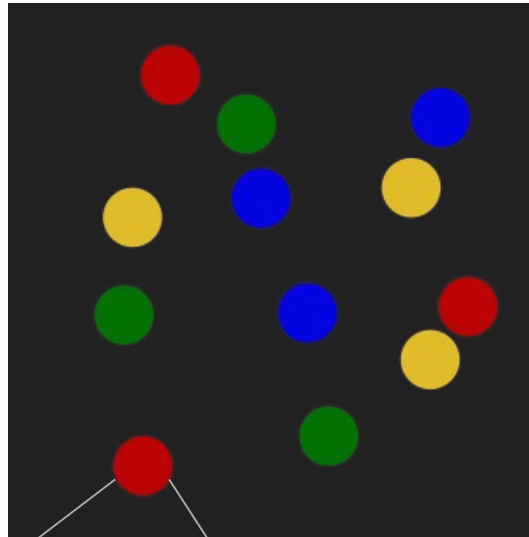

in phase

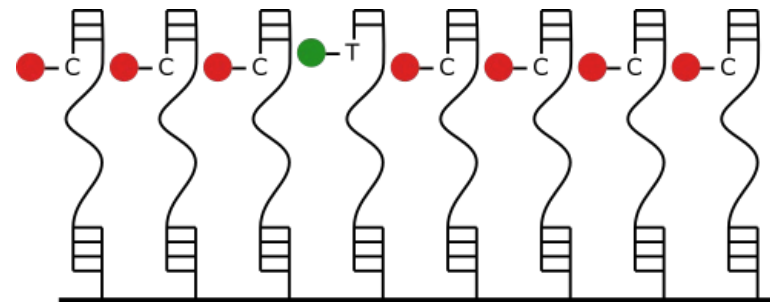


out of phase

Long reads (3rd generation)


- Crux: maximizing signal from a single-molecule base read-out
- Single molecule, so no out-of-phase signal
- Two frequently used platforms:
 - PacBio SMRT sequencing
 - Oxford Nanopore Technology

Oxford Nanopore technology



- Based on changes in electrical current
- Well-known for its scalability and portability
- ~95-97% accuracy

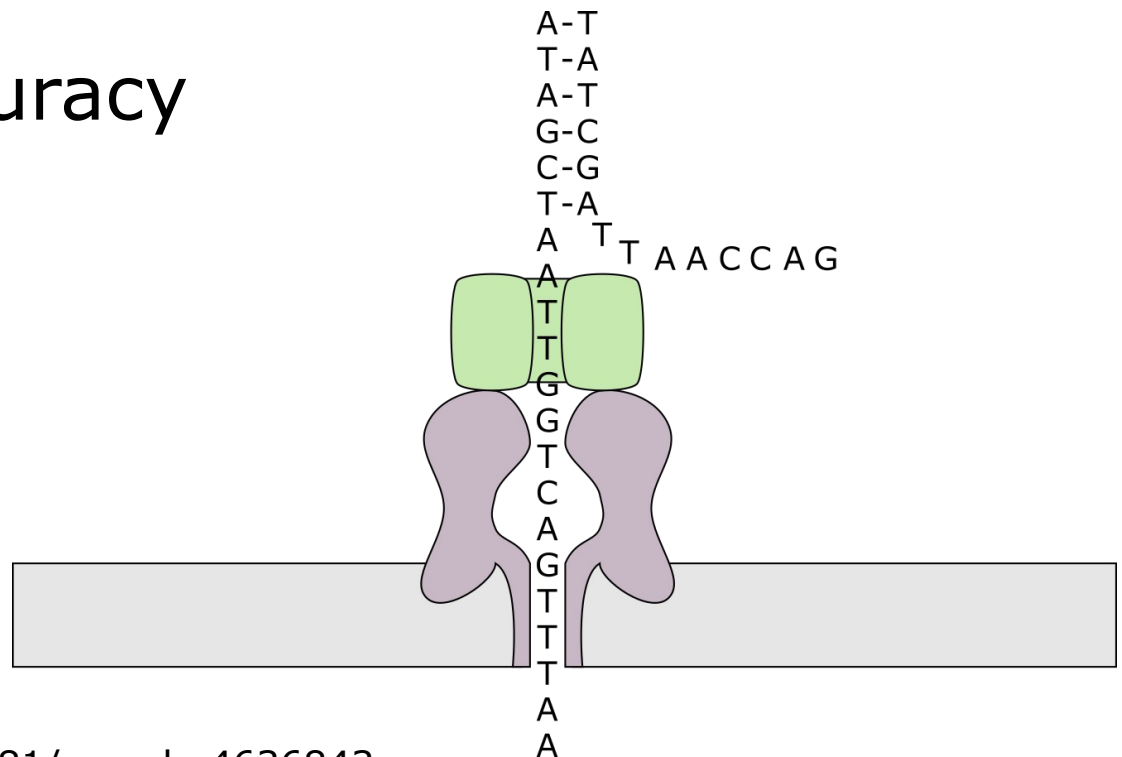

PacBio sequencing

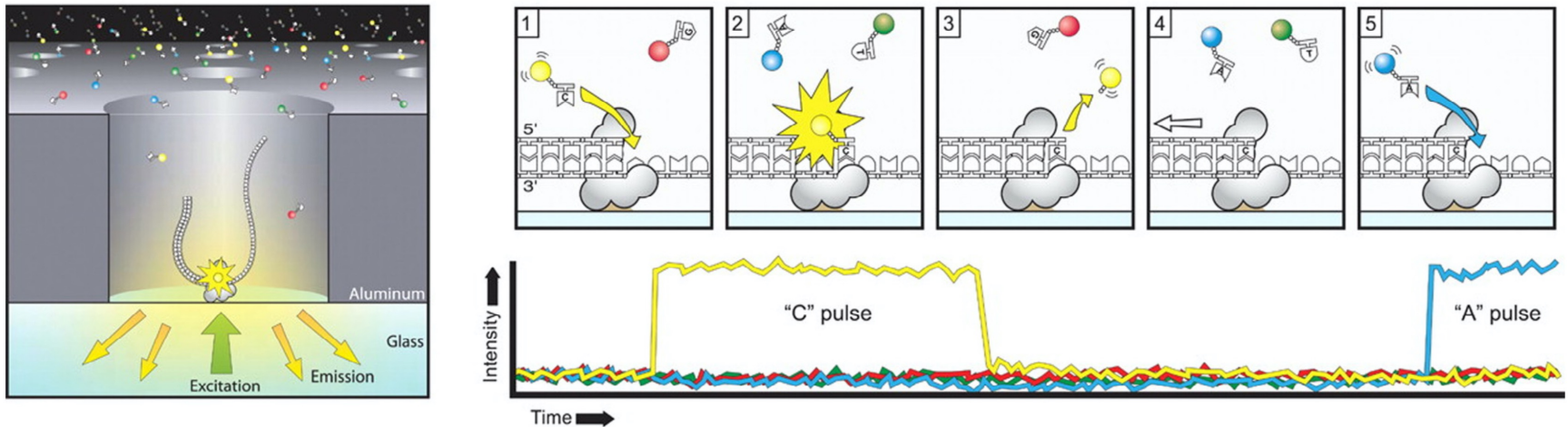

- Polymerase bound to ZMW bottom
- Circular molecules
- Single read out $\sim 90\%$ accuracy
- CCS (HiFi): single molecule sequenced multiple times

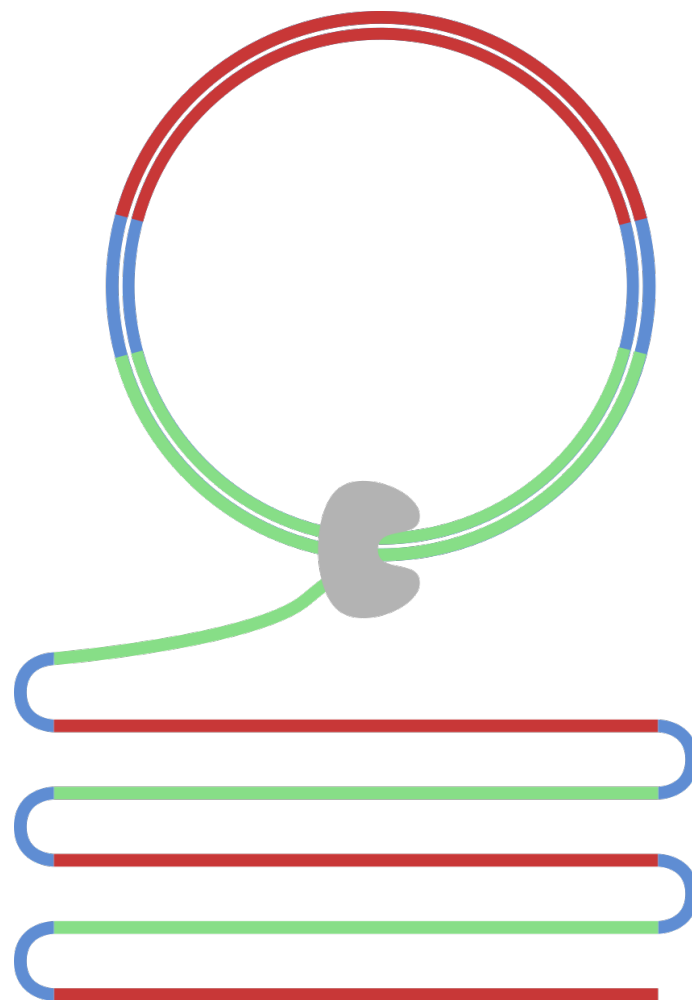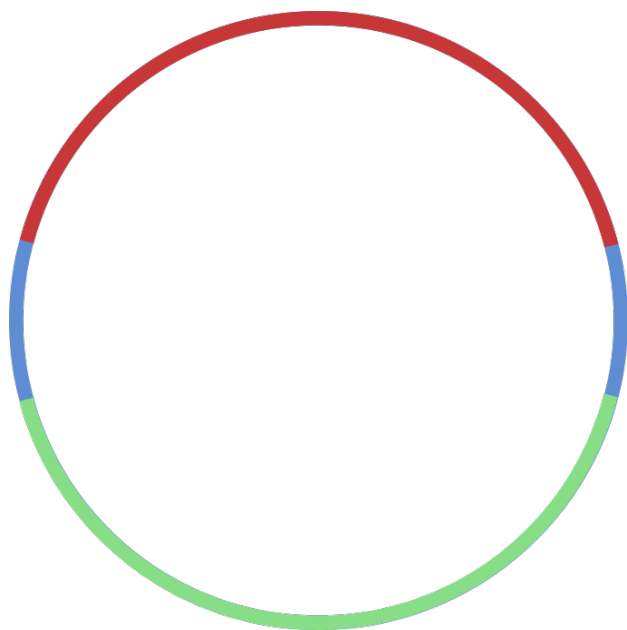

Hi-Fi read

Quiz Question 6 and 7